

ZELOSTAT[®]

Needle[™]

CE 0344



ZELOSTAT[®]

A new experience for injection

Super thin-walled needles ensure smooth and comfortable fluid flow while specially designed needle hub with ultra low dead space expels residual fluid to a minimum. Once you try ZELOSTAT Needles, the smooth needle insertion and fluid flow will be an unforgettable experience and you won't be able to go back.


SUPER THIN WALL

Allowing smooth and comfortable injection. For patients, injection with ZELOSTAT needles conduce comfortable treatment experience.

Comparison of Diameters


Nominal Gauge	Product	OD (mm)	ID (mm)
31G	 ZELOSTAT	0.25	0.20
	Conventional Needles	0.26	0.15
34G	 ZELOSTAT	0.18	0.12
	Conventional Needles	0.18 or 0.20	0.10
35G	 ZELOSTAT	0.16	0.105
	Conventional Needles	Not available	Not available

* Reference value

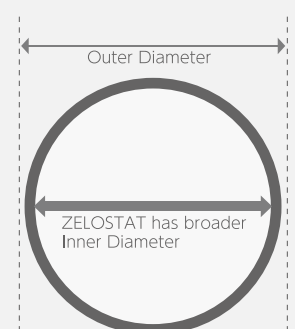


ULTRA LOW DEAD SPACE

Specially designed needle hub with ultra low dead space expels residual fluid to a minimum. ZELOSTAT needles minimize injectable waste

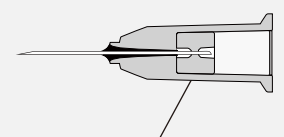

Elastic body infills dead space



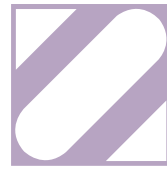

ZELOSTAT

Needle™

35G IS IN THE LINE

ZELOSTAT 35G needles enable an injection with less pain, less bleeding and less invasive you never experienced. While we realize the outer diameter of the needle tube as 0.16mm, its inner diameter remains 0.105mm which is even broader than conventional 34G, typically 0.10 mm.



APPLICATIONS


Size	Use
25G×25mm	Withdrawing injectables
31G× 4mm	Mesotherapy, PRP, B. Toxin, HA (low-middle viscosity)
31G× 9mm	
31G×13mm	B.Toxin, HA (low-middle viscosity)
34G× 4mm	Mesotherapy, PRP, B.Toxin
34G× 9mm	
35G× 4mm	Mesotherapy, PRP, B. Toxin (under eye, lip)

made in Japan

Item Code	Gauge (Outer Diameter)	Length	Package
NZ-25S-25-100	25G (0.50mm)	25mm	100pcs
NZ-31S-04-100	31G (0.25mm)	4mm	
NZ-31S-09-100	31G (0.25mm)	9mm	
NZ-31S-13-100	31G (0.25mm)	13mm	
NZ-34S-04-100	34G (0.18mm)	4mm	
NZ-34S-09-100	34G (0.18mm)	9mm	
NZ-35S-04-100	35G (0.16mm)	4mm	50pcs
NZ-35S-04-050	35G (0.16mm)	4mm	

Needle Hub Colors According to ISO 6009:2016

For use exclusively with luer lock syringes.

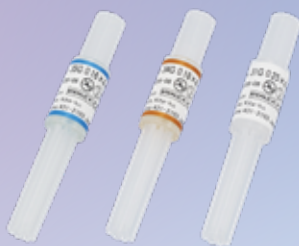

Plastic Package



Needle samples are available.

www.asti-med.com

Mail: medical@asti.co.jp